



Go with the flow...

Evolve Flow bondable tubes were designed in conjunction with leading orthodontists to incorporate all of the features needed for enhanced performance. Available in Roth* and MBT* .018 and .022.

Low Profile Hook

Low profile hooks guarantee patient comfort and allow for easy placement of attachments.

Identification Marks


Each tube is colour coded on the buccal surface for easy molar quadrant identification.

Distal Notch

Evolve Flow features a distal notch for ligation and tieback of auxiliaries.

Low Profile

The contoured surfaces are low profile and tapered for patient comfort.

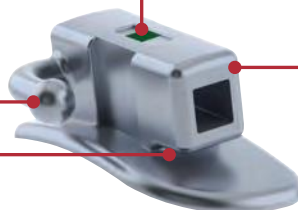

Contoured Base

Anatomically contoured base with occlusal indent ensures accurate fit, whilst the 80 gauge mesh pad ensures superior bond strength.

The Funnelled Entrance

The Evolve Flow funnelled entrance allows for easy wire insertion.




Positioning Guide

Occlusal / gingival positioning guides ensure a secure grip and ease of placement.


MBT* .022	1st Molars	2nd Molars
Upper Left	DB07-1360UL	DB07-1372UL
Upper Right	DB07-1360UR	DB07-1372UR
Lower Left	DB07-1360LL	DB07-1372LL
Lower Right	DB07-1360LR	DB07-1372LR

Roth* .018	1st Molars	2nd Molars
Upper Left	DB07-1370UL	DB07-1373UL
Upper Right	DB07-1370UR	DB07-1373UR
Lower Left	DB07-1370LL	DB07-1373LL
Lower Right	DB07-1370LR	DB07-1373LR

1st Molar Evolve Flow Tube

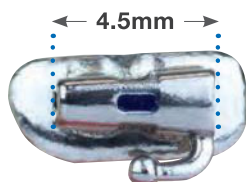

2nd Molar Evolve Flow Tube

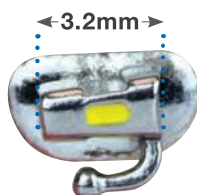

Roth* .022	1st Molars	2nd Molars
Upper Left	DB07-1350UL	DB07-1371UL
Upper Right	DB07-1350UR	DB07-1371UR
Lower Left	DB07-1350LL	DB07-1371LL
Lower Right	DB07-1350LR	DB07-1371LR

*The prescriptions are not implied to be an exact version of any other system, nor do we claim any endorsement by the doctor(s).